


RT-SPOT[®]


Beekley
Medical[®]
when your diagnosis
must be right[®]

The only skin marker for CT simulation that is elevated from the skin line and helps to exclude the marker from the automatic body contour


- ▶ Minimizes dose perturbation and helps ensure dose accuracy
- ▶ Saves time by reducing the need for manual contouring and density override
- ▶ Flexible and easier to form into desired shapes
- ▶ Stays in place on the skin and does not lift
- ▶ Lower density, non-metallic line, minimizes artifact

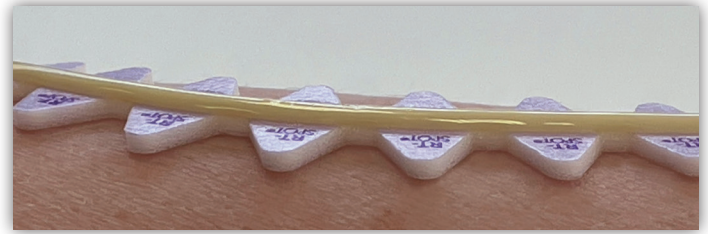

▶ RT-SPOT's medical grade foam backing elevates it from the skin line

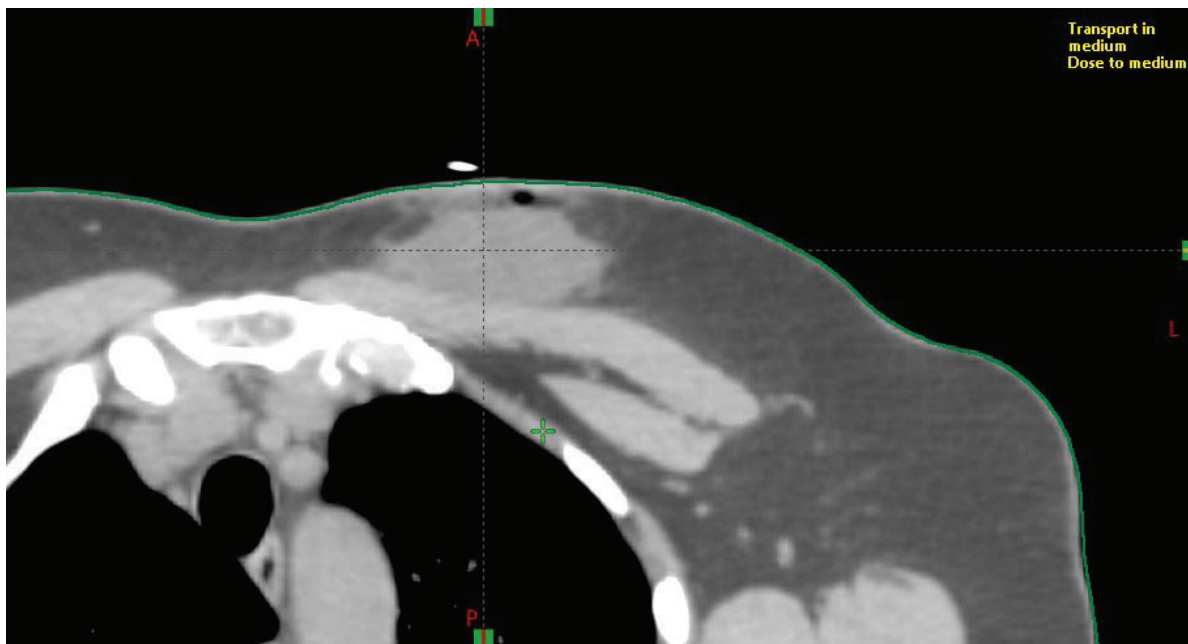

◀ Axial CT Scan of female breast radiation patient with surgical scar identified by RT-SPOT. RT-SPOT excluded from body contour.

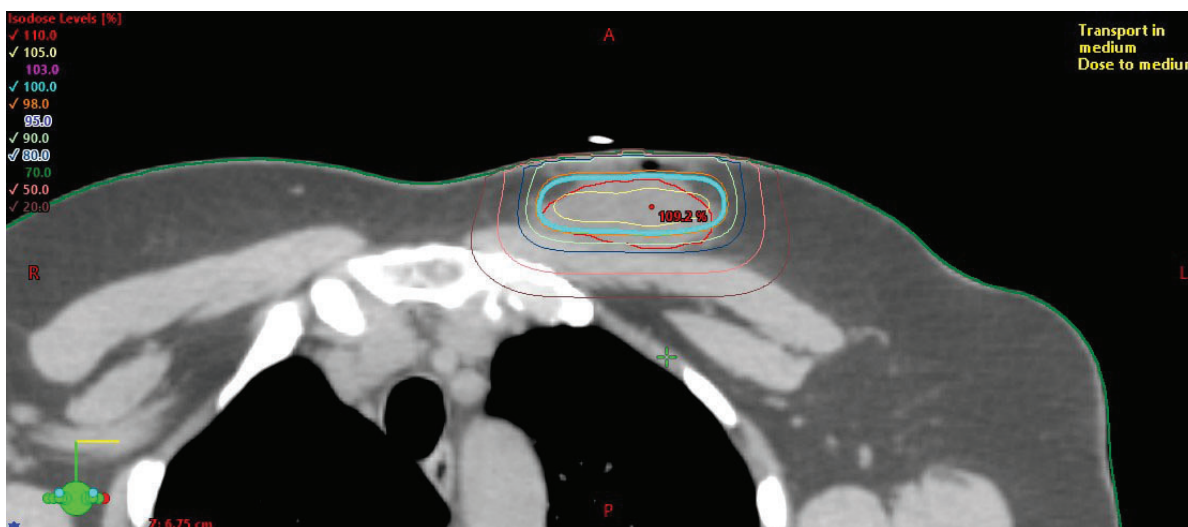

◀ Corresponding image of electron dose distribution unperturbed by RT-SPOT.

RT-SPOT[®]

Beekley Medical[®]
when your diagnosis
must be right[®]

Existing skin markers affect dose calculations and must be manually contoured out, increasing the dose planning time per patient.


Manual contouring can be imperfect leading to dose inaccuracies or rework of treatment plans.

RT-SPOT is an elevated skin marker designed to minimize dose perturbation and reduce the need for manual contouring.

RT-SPOT is a flexible line that contours more easily into desired shapes and stays in place without lifting. Ideal for marking field borders, tangents, scars, sarcomas and larger treatment areas.



▲ The moldable backing allows RT-SPOT to bend into various shapes – even small circles

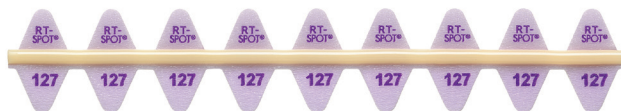


REF 127 **QTY** 140cm cut to measure roll

- 2.0mm diameter, elevated, non-metallic line
- For marking field borders, tangents, and scars

1.800.233.5539 • info@beekley.com • beekley.com



Customers outside the U.S. – Contact your local distributor for pricing and product availability.

To locate a distributor, call +1.860.583.4700 or email international@beekley.com

 Manufactured by Beekley Corporation
One Prestige Lane, Bristol, CT 06010 USA
Tel: 1.800.233.5539 or +1.860.583.4700
Fax: 1.800.735.1234 beekley.com

BEEKLEY, BEEKLEY MEDICAL, SPOT, and WHEN YOUR DIAGNOSIS MUST BE RIGHT are Reg. U.S. Pat. & Tm. Off. BEEKLEY and SPOT are Registered Community Trademarks. BEEKLEY is a registered trademark in Canada. RT-SPOT is an unregistered trademark of Beekley Corporation.
© 2021 Beekley Corporation. All rights reserved.

WARNINGS AND PRECAUTIONS:

- Do Not Use in Skin Folds
- Do Not Use Under Physical Bolus
- Single Patient Use Only; Do Not Reuse 
- External Use Only; Do Not Ingest or Swallow
- Do Not Place on Compromised Skin
- Remove After Use
- Not Made With Natural Rubber Latex
- Lead-Free
- Non-Sterile 

If any serious incident occurs in relation to this device, report it to Beekley Medical and to the competent authority of the Member State in which the user and/or patient is established.



EMERGO EUROPE
Prinsessegracht 20
2514 AP The Hague
The Netherlands



REV: RTS_INF_0821